

# **Visitors to Care Homes in Northern Ireland**

## **Easy Read Guide:**


# **COVID-19 self-testing using Lateral Flow Devices (LFDs)**

## **How to use a home-testing kit**

# About this booklet



**Coronavirus** or **COVID-19** is an illness affecting people all over the world.

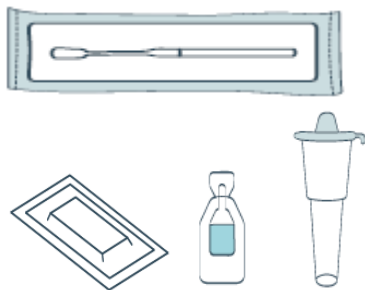

You can use a **home testing kit** to find out if you have COVID19. This booklet tells you how to use the kit.



A **home testing kit** is a pack with everything you need to do the test yourself at home. You can ask a friend/relative to help you with the test if you need to.



Some people have COVID19 but don't know it. By using these tests we can find out who has COVID19 and stop it spreading.

# About the test



The test will take about 15 minutes. The test result will be ready 30 minutes after you finish the test.

The result means the answer to the test, and tells you if you have COVID-19.



You need to

- A. Get ready
- B. Get the test ready
- C. Take a sample with a swab
- D. Test the sample
- E. Get your result
- F. Tell us your result
- G. Throw your test kit away



This booklet tells you how to do the test. You must read each part of this booklet and do what it says.

# Before you start the test



Keep the kit at normal room temperature. This means it should not be left in the sun, or in a fridge or freezer. Keep the kit away from children.



Use 1 kit for each person. You can only use each part of a kit once. Don't use the kit on animals.



If your nose is pierced on both sides, take the piercing out on 1 side. You will need to take a sample from the side with no piercing.



Don't eat or drink for 30 minutes before you use the kit. If you had a nosebleed in the last day, wait 24 hours before you use the kit.

## A. Get ready

1



Read this booklet carefully.  
You can also watch a video  
showing how to take a test by  
clicking on this link

[www.gov.uk/covid19-self-test-help](https://www.gov.uk/covid19-self-test-help)

2




Get a flat area ready, like a  
table. Make sure it is clear,  
clean and dry.

3



Wash your hands for 20  
seconds. Use soap and warm  
water, or hand sanitiser. Dry  
your hands.

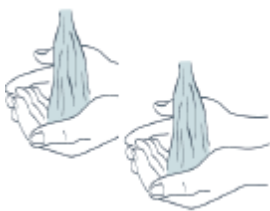

If you are doing more than 1  
test, clean the table and wash  
your hands between each test.

---

4



Make sure you have everything you need. If anything is missing or broken do not use the kit.

### From the kit



Test strip in a closed pack




Swab in a closed wrapper



Extraction buffer sachet



Extraction tube

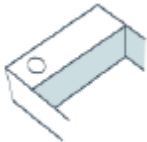

Extraction tube holder (part of the box)



Waste bag

### From your home



a watch or clock



a mirror




tissues



hand sanitiser or soap and warm water

## B. Get the test ready

5

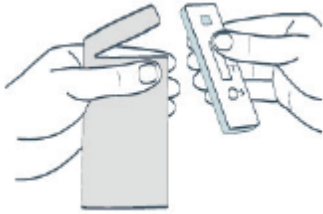


Take the test strip out of the closed pack. Put the test strip on your table, and the pack in the bin. You must start the test in the next 30 minutes.

6




Open the sachet carefully by twisting it or snapping it. Be careful not to spill any of the liquid in the sachet. Keep it away from your face.

7

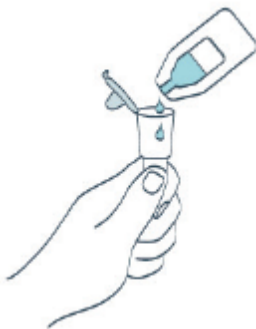

Open the extraction tube. Slowly squeeze all of the liquid in the sachet into the tube. Try not to let the sachet touch the tube.

8

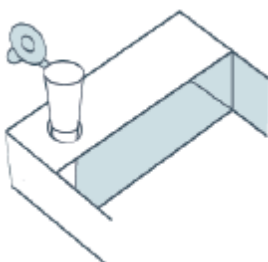

Put the extraction tube in the tube holder. This is a hole inside the box, it makes the tube stand up so you don't spill the liquid inside it.

9




Blow your nose into a tissue and throw the tissue in the bin.

10



Wash and dry your hands again. If you are helping someone else do a test, you must both wash your hands.

11

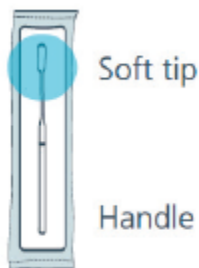

Pick up the swab in the closed wrapper. Make sure you know which end is the soft tip of the swab.

12



Open the closed wrapper and take out the swab by holding the handle. Never touch the soft tip with your hands.

## C. Take a sample with a swab

13



Hold the handle of the swab with your fingers. Open your mouth wide. Rub the soft tip of the swab over both tonsils 4 times.

**Tonsils** are soft pads at the back of your throat, 1 on each side. If you have no tonsils rub the swab on the sides of your throat.

Do not let the swab touch anything else, like your tongue or teeth.

14




Take the swab out of your mouth and put it inside 1 side of your nose. Push it up until it stops. Rub the swab in a circle around the inside of your nose 10 times.

When you use the swab it might feel strange but it should not hurt.


## D. Test the sample

15

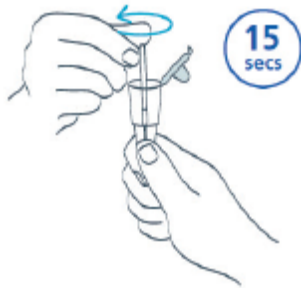


Pick up the **extraction tube** and open the lid. Put the soft tip of the swab inside the tube in the liquid. Stir the liquid with the swab and wipe the swab against the side of the tube for 15 seconds.

16

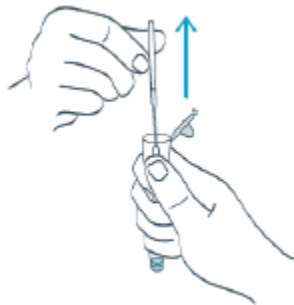

Take the swab out of the extraction tube. When you do this squeeze the tube onto the swab to help get all the liquid out of the swab. Put the swab in the **waste bag**.

17

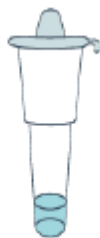

Put the lid back onto the extraction tube. Make sure the lid is on tight so that it cannot leak.

18



Turn the extraction tube upside down and squeeze 2 drops of the liquid onto the **test strip**. Make sure the drops are in the **specimen well** which is the hole next to the letter S.

19



Look at your clock or watch.  
You need to wait for 30  
minutes for your test to be  
ready. Work out what time your  
test will be ready.



Don't move the test strip. It  
must stay flat and still.



While you wait, some lines  
might show on the test strip.  
Keep waiting as the lines might  
change.




After 30 minutes your test will  
be ready and you can read the  
result.




On the next page we will tell  
you what the lines mean and  
what your test result is.

## E. Get your result

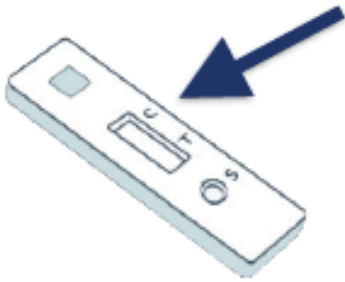

Look at your test strip and see if there are any lines next to the letter C or the letter T. There are 3 possible results. Match yours to a picture on this page to find your test result.

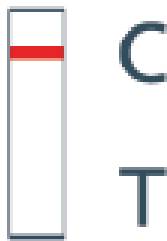

### Negative result

If you have 1 line by C, and no line by T this is called a negative result. You probably do not have COVID-19. You can do another test if you get ill.



### Positive result

If you have 1 line by C, and 1 line by T this is called a positive result. The lines can be bright or faint. You have COVID-19.



### Void result

If you have no lines, or 1 line by T this is called a void result. This means that the test did not work. You need to do another test.

## F. Tell us your result



You must tell us the result of your test. This will help us learn more about COVID-19 and help more people.



To tell us your test result you will need the QR code or the number on your test strip. The **QR code** is a pattern of dots that a computer or mobile can read.



You can tell us your result by going to this website [www.gov.uk/report-covid19-result](https://www.gov.uk/report-covid19-result)



Or you can tell us your result by phoning us between 7am and 11pm.

In England, Wales and Northern Ireland call this number **119**

In Scotland call this number **0300 303 2713**

## G. Throw your test kit away



After you have done the test, put all parts of the kit in the waste bag.



Put the waste bag in your bin.

# What next?



If you have COVID-19 you and everyone who lives with you must stay at home. You must follow the rules in the area where you live.



If you are not better after 7 days call 111. If you are very ill call 999.



If you need more help with the test visit this website [www.gov.uk/covid19-self-test-help](https://www.gov.uk/covid19-self-test-help)

The full version of this document is available at the link below and is called

**“Visitors to Care Homes in Northern Ireland - Your step-by-step guide for COVID-19 self-testing using Lateral Flow Devices (LFDs)”**

[Testing for visitors to care homes.](#)

Please read the full document for more information